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EN | USAGE INSTRUCTIONS

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INSTRUCCIONES DE USO DE
IMPLANTES
IMPLANT USAGE
INSTRUCTIONS
INSTRUCTIONS D'USAGE
DES IMPLANTS

ZES-CA-D-1146-G-Rev-04
Noviembre 2020

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KYT®

Implantes dentales KYT®

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INFORMACIÓN IMPORTANTE.

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USAGE INSTRUCTIONS

IMPORTANT INFORMATION.

Read this prospectus carefully before using the product.

This document contains basic information for the use of KYT® dental implant systems. This documentation has been drafted as a quick guide to be consulted by the final user of the system and does not provide the indications and technical specifications necessary to properly use the products in the KYT® dental implant systems. This is not an alternative nor substitute for specialized training or professional clinical experience.

The KYT® dental implant systems must be used following a strict surgical and prosthetic protocol that does not excessively deviate from any other conventional implant on the market. As with these systems, adequate treatment planning is necessary. However, one must take some differences from other implant systems into account, both in planning treatment and in the surgical and prosthetic protocol in order to avoid errors in the working process and manipulation of the products.

Read carefully the drilling sequences and the general system recommendations specific to each product, the instructions for use and maintenance before using a product of the KYT® implant system. These can be found on our website www.kyimplantesdentales.com or requested from your nearest KYT® distributor.

1. Description of the system.

KYT® has developed several dental implant systems in order to provide dental treatments in accordance with the existing scientific literature and current clinical standards.

KYT® develops and markets ultimate generation surgical instruments and a wide array of restoration solutions to meet each patient's particular needs.

2. Description of the products.

Implants

KYT® dental implants are made of Grade 4 commercially pure biocompatible titanium in accordance with ASTM/F-67 and ISO 5832-2, with a specially treated surface to improve osseointegration. The dental implants may have different diameters, lengths, shapes and connections depending on their design.

3. Usage indications.

Before performing any treatment with KYT® dental implants, it is essential that you evaluate the patient and the patient's possible contraindications for treatment with dental implants (clinical diagnosis, radiology and study of dental models). Proper treatment with dental implants must first have a proper prosthodontic rehabilitation design that meets the patient's needs and expectations.

Implants

KYT® implants are indicated for both the upper and lower jaw and are designed to restore normal chewing function and dental aesthetics. They come in various lengths and diameters to suit different clinical situations.

Important. Carefully read the clinical indications of each implant in the catalogs and surgical protocols of each product before carrying out surgical planning. Do not use KYT® dental implants in sites or situations which are not indicated or described as "Indicated" in their catalogs or surgical protocols.

4. Contraindications.

A preoperative medical examination of the patient is necessary to determine if there are contraindications or risk factors for treatment with dental implants. KYT® dental implants should not be used in patients who lack the minimal medical conditions necessary to carry out prosthodontic treatment and rehabilitation with dental implants.

Relative. The clinical professional is responsible for evaluating the potential benefits and risks for the particular treatment: Bone deficiency - Bone quantity/quality- Patient's age- Patient's health- Poor oral hygiene- Localised infections at implant site- Residual roots at implant site- Pregnancy- Stress- Mental factors- Prosthetic heart valves- Bruxism- Peridental pathology- Gum pathology- Tobacco use- Alcohol abuse- Drug Addiction- Addiction in general, among others.

Absolute. The manufacturer does not recommend the treatment in the next cases: Cardiovascular diseases or disorders - Haematological diseases - Maxillary sinusitis- End-stage disease- Osteoporosis or bone disorders- Epilepsy - Neoplasms treated with chemotherapy - Blood dyscrasias - Cancer - Severe or poorly controlled diabetes- Untreated periodontitis - Infectious diseases - Leukaemia - Immune system alterations- Allergy to titanium- Allergy stainless steel. - Radiation therapy treatment- Corticosteroid treatments - Anticoagulant treatments - Involuntary muscle disorders - Severe liver disorders - Endocrine disorders- Uncontrolled hypertension- Psychiatric and/or psychological disorders.

5. Storage.

KYT® products must be stored in their original packaging at a temperature between -10 and 55°C in a dry, clean location protected from sunlight and adverse conditions.

6. Packaging and manner supplied.

KYT® dental implants are presented in a cardboard box protected by a transparent shrink-wrapped plastic and with a label identifying the product, containing:

- These explanatory Instructions for Use of the product.

- Heat-sealed blister with a sticker indicating the manufacturer, the complete product reference and the batch number. Inside the package is the plastic container within the implant carrier that holds the implant and its transfer, the final screw and the closure screw of the implant.

- Stickers with the necessary information for the correct traceability of the implant. One will be attached to the Patient's Implant Card and another will be kept in the clinic's register.

KYT® guarantees that all of its products follow a strict manufacturing, quality control and cleaning process before being packaged. Before using a KYT® implants products, you must check for proper integrity of the packaging and ensure that it is not damaged. If there is a defect or damage, the product must not be used, must be discarded and a replacement product ordered. KYT® will replace any damaged product at no cost so long as it does not show signs of previous handling.

7. Cleaning and sterilization.

KYT® dental implants are sterilised using 25 kGy beta-radiation. If the sterile container is altered, mistreated, damaged or broken, sterility is not guaranteed and the product must not be used. Do not use the product after the expiry date shown on the product container.

Resterilization of KYT® dental implants is NOT permitted.

Important. The heat-sealed blister and the bottle containing the container with the implant are NOT sterile on the outside, only the inside of the bottle is sterile; do not place the blister or the bottle in the sterile surgical field. This sterile container must NOT be opened until just before placement of the implant.

8. Precautions.

Surgery and Instruments

The drilling sequences and general recommendations of the system describe in detail the precautions and important factors to be taken into account during the surgical process of placing the dental implants. Prior to the operation, it is essential to have consulted the drilling sequences and general system recommendations of the product, as well as all the indications on its correct use indicated.

Due to the average dimensions of the products used, among the most important measures you must take special care to ensure that these products are not accidentally ingested or aspirated by the patient. All KYT® instruments for manual use incorporate retention elements for use with dental thread or silk to avoid accidental ingestion.

You must do everything possible to minimize the surgical time and damage to the area where the implants will be placed, paying special attention to bone trauma and asepsis of the surgical field. Proper preparation of the surgical bed is a fundamental factor for successful treatment. This process requires that you follow the surgical burring protocol provided by the manufacturer closely for each type of implant and proper use of the surgical burrs, both in the number of uses and the working speeds. The surgical burrs must be used with constant abundant irrigation to avoid bone overheating.

Pay special attention to not exerting excessive forces or torque during placement of the implant. It is important to control the insertion torque of the implant, both in mechanical insertion and in manual insertion.

Attention. Use of the surgical burrs, Maximum 45 Uses. Surpassing the number of uses of the surgical burrs puts the osseointegration process and therefore the dental treatment at risk.

The maximum insertion torque for the dental implants is 50 Ncm. Surpassing the maximum insertion torque for the implants may cause severe damage in the dental implant, its connection, the carrier and the clinical screw.

9. Warnings.

Treatment planning and placement of KYT® dental implants require specific dental training. The product specifications themselves are not guarantee of proper use. We recommend that clinical users take courses with specific theoretical/practical training in order to understand the proper working techniques and protocols for each product, including the biomechanical, radiographic and prosthodontic requirements associated with implant treatment. The end user must be knowledgeable and inform himself about the current state of the technique for each product and its possible applications.

Before using KYT® implant systems or any of their prosthodontic attachments or components, it is necessary to understand and be familiar with the corresponding surgical and prosthodontic procedures. The patient must meet some sufficient anatomical and psychological conditions to undergo treatment with dental implants. The person responsible for implant treatment must

watch for proper treatment planning, guaranteeing the safety of the treatment with a minimal margin of error, respecting the patient's vital mouth structures and general health.

The procedure for the use of surgical drills, taps and other surgical instruments necessary for the placement of the implants is specified in the drilling sequences corresponding to each implant. Placement of the implant and prosthodontic planning must adapt to the patient's individual conditions, particularly to proper distribution of forces. Passive adjustment in the prosthodontic rehabilitation, adjustment of occlusal contact points at the opposing maxillary and avoiding the onset of excessive lateral forces must be achieved. Treatment with an insufficient number of implants, improper choice of the size or an inappropriate position to support and transmit the expected loads, may lead to mechanical failure of the implant, the abutment or the prosthetic attachment, the clinical screws due to overload or fatigue and substantial loss of bone around the implant.

Lack of an adequate amount of residual bone, onset of infection or diseases in general and changes in the patient's habits are some potential causes of osseointegration and treatment failure. Lack of bone or soft tissue may lead to unfavorable insertion of the implant and lead to a subsequent deficient aesthetic outcome. Inadequate prosthodontic rehabilitation may cause failure of the entire rehabilitation treatment.

Reuse of single-use products carries the risk of possible deterioration of their characteristics. The geometry of the product may have suffered damage from previous use or due to inadequate use. This involves the risk of surgical or prosthodontic failure and/or damage to the patient's health.

Electrosurgery is not indicated due to the conductivity of dental implants.

10. Adverse or side effects.

Current implant technology can have several adverse or undesired effects that are documented in the specialized scientific literature published in journals or specialized books on dentistry.

However, the most relevant are:

Postoperative malaise- Local pain- Local inflammation- Hematomas- Local or systemic infections- Difficulty speaking- Bone loss and/or fractures- Loss of the implant- Damage to adjacent teeth- Fractures of implants and prosthodontic components- Instrument fracture- Damage to the tooth nerve- Damage to the maxillary sinuses, among others.

11. Information on responsibility, safety and guarantee.

Our usage and application indications for our products are based on the internationally published literature, current clinical standards and our clinical experience with our dental implant system, so they must be understood as general information on indications. Handling and use of the KYT® dental implant systems, being out of the control of KYT®, are the sole responsibility of the end user of the products. KYT Implantas Dentales, S.L. and the distributor of the KYT® implant systems decline all civil responsibility for possible damage or injury caused by poor handling of the product or any other situation not considered in the protocol for proper use of its products.

The end user must ensure that the product used is adequate for the planned purposes and procedures. Neither these usage instructions, the working protocols or manipulation of the products release the user from this obligation. The KYT® dental implant systems must be used, handled and applied by professionals with the necessary qualifications required according to current legislation of each country. Total or partial use, handling or application of any phase of the KYT® dental implant systems by non-qualified personnel or by personnel who do not have the necessary training will automatically void any type of warranty and may cause serious damage to the patient's health.

The KYT® dental implant systems products are part of their own system, with their own design characteristics and working protocols that include dental implants, prosthodontic attachments or components and surgical or prosthetic instruments. The use of original KYT® products in combination with inappropriate elements or components or products from other manufacturers may lead to failure of the rehabilitation treatment, cause tissue damage, lead to damage to bone structures, inadequate aesthetic outcomes and severe damage to the patient's health.

For this reason, ONLY original KYT® products should be used. The clinician is responsible for ensuring that original KYT® products are used and that they are used in accordance with the corresponding instructions for use and drilling sequences throughout the entire process of implant treatment. The use of components, instruments or any other non-original KYT® product, whether used singly or in combination with any of the original products, will automatically result in the loss of any guarantee, either total or partial, on the KYT® products.

To make use of the guarantee provided for original KYT® products, all products used in implant treatment must be original KYT® products and the requirements for the intended use of the product, the original KYT® instructions for use and the original KYT® drilling sequences must have been complied with. Both the person responsible for the dental treatment and the patient must provide all information required by KYT® for study. All original products that have been used to treat patients with risk factors or contraindications expressed in this document are excluded from the total or partial guarantee.

Warning. Not all KYT® products are for sale in every country.

For further information, please consult the specific product catalogues on our website www.kytimplantesdentales.com or request them from your nearest KYT® dealer.

KYT® and all their logos are registered trademarks of KYT.

The KYT® implant systems are a registered trademark, as are some of the products mentioned or not mentioned in this catalogue. KYT Implantés Dentales, SL reserves the right to modify, change, eliminate or update any of the products, prices or technical specifications referenced in this document or in any of its catalogues without prior notification. All rights reserved. Any complete or partial reproduction of this documentation, in any medium or format, is forbidden without the corresponding written authorization from KYT Implantés Dentales, S.L.

12. Magnetic resonance.

The KYT® dental implant systems products have not been evaluated for safety and compatibility nor have they been tested for heating and migration in the magnetic resonance environment.

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Rx Only

Caution: federal law prohibite dispensing without prescription.

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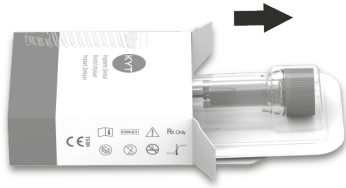
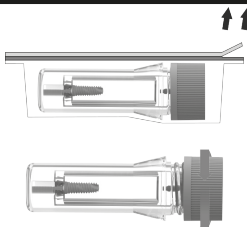
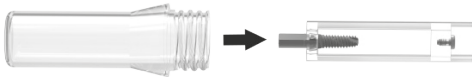
EN Carefully remove the seal from the implant blister.

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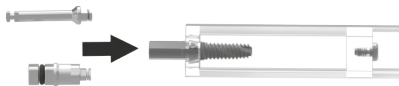
EN Drop the vial carrying the implant onto a sterile cloth in the surgical area.

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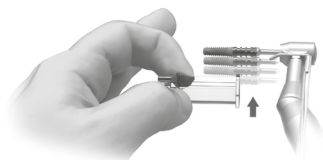
EN Hold the vial carrying the implant with one hand and insert the mechanical or manual adapter into the transfer.

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EN Remove the implant transfer unit by pulling upwards vertically towards the vial.

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